

Invacare® Perfecto₂™ Oxygen Concentrator




Perfecto₂™



Yes, you can.®



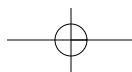


Invacare® Perfecto₂™ Oxygen Concentrator

The new Perfecto₂™ is the smallest, lightest, quietest and most efficient 5-liter concentrator ever produced by Invacare. The high performance Perfecto₂ unit incorporates all the best features of the market-leading Platinum™ XL concentrator, while significantly improving upon the key specifications of the Platinum XL and offering a contemporary, patient preferred design and color. Available with or without the SensO₂® oxygen sensor, the Perfecto₂ concentrator represents a quantum leap forward in the evolution of oxygen concentrators.



Yes, you can.®

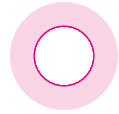

Invacare® Perfecto₂™ Oxygen Concentrator

Model no. IRC5P 5-liter concentrator
 Model no. IRC5P02 5-liter concentrator w/ SensO₂® (shown below)

NEW Features

- 75% quieter than Platinum XL
- 25% more energy efficient than Platinum XL
- 17% lighter than Platinum XL
- 33% smaller than Platinum XL
- Contemporary design and color

Key Features

Top handle	For ease of mobility
Servicing parts compatibility	Compatible with many Platinum concentrator parts
Alarms and safety systems	Battery-free power-loss alarm; Sieve-GARD™ monitors sieve performance; SensO ₂ ® oxygen monitor; Compressor 35 psi pressure-relief valve, low flow alarm
Pressure-based system	Superior to timing-based systems
Self-diagnostic electronics	Minimizes troubleshooting time
Simplified P.M. schedule	Replace P.M. components as needed

* Based on data gathered from independent, third party repair centers.

Specifications

Dimensions	
Width	13"
Height	23"
Depth	11.5"
Weight	43 lb.
Shipping weight	47 lb.
Sound level	43 dB
Average power consumption	280 watts @ 3 LPM
Number per pallet	27 units
Oxygen output concentration levels	95.6% to 87% at all flow rates
Maximum outlet pressure	5 +/- 0.5 psi.
Flow range	1/2 - 5 L/min.
Electrical requirements	120 VAC +/- 10%, 60Hz
Rated current input	4.0 A
Pressure relief mechanism operational	35 psi +/- 3.5 psi
Filters	Cabinet, outlet HEPA, compressor inlet
Limited warranty	5 years



Invacare® Perfecto₂™ Oxygen Concentrator



HomeFill™ connection.



Filter access and adapter storage.

The HomeFill™ compressor and Ready-Rack also have been updated to coordinate with the Perfecto₂.

Model no. IOH200-BK HomeFill compressor
 Model no. IOH270 Ready-Rack



Model Selection Guide

Model No.	Units/Case	Limited Warranty	Description
IRC5P	1	5 Years	Perfecto ₂ 5-liter concentrator
IRC5P02	1	5 Years	Perfecto ₂ 5-liter concentrator with SensO ₂
IOH200-BK	1	3 Years	HomeFill Oxygen Compressor (Graphite color)
IOH270	1	1 Year	HomeFill Ready-Rack (for use with Perfecto ₂ only)

Invacare Corporation
 www.invacare.com

USA

One Invacare Way
 Elyria, Ohio
 44035-4190
 (800) 333-6900

Canada

570 Matheson Blvd. E.
 Unit 8
 Mississauga, Ontario
 L4Z 4G4 Canada
 (800) 668-5324

© 2007 Invacare Corporation.

All rights reserved.
 Trademarks are identified by the symbols ™ and ©.
 All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.

Specifications are subject to change without notification.

©2007
 Invacare Corporation
 Form No. 07-007
 Rev. 05-07

