



Lightest Foam/Gel Overlay

- A new standard in pressure management
- Extremely light-weight for easy handling
- Mattress overlay incorporates a Solid Gel patterned layer that is inherently visco-elastic. It works in responding to pressure by deforming slightly, thereby dispersing critical weight over a large surface area
- Non migrating gel is engineered into a high density 2" solid foam base with a specially designed 1" convoluted top to provide maximum comfort and pressure management
- Overlay reduces body interface pressure and provides a comfort at critical points over a wide range of body shapes and sizes
- Since the gel will not migrate or pool, the overlay is effective in any bed position
- Our non-aqueous gel will not migrate or leak like water-based gel products
- Overlay delivers same consistent results over the lifetime of the product
- Light weight feature reduces shipping costs
- Cover is fluid resistant, anti-bacterial and anti-fungal. Meets CAL 117
- 2 year limited warranty


103680	36"x 80"x 3", 26 lbs, 275 lbs capacity	1 ea
103684	36"x 84"x 3", 26 lbs, 275 lbs capacity	1 ea
104280	42"x 80"x 4", 34 lbs, 400 lbs capacity	1 ea

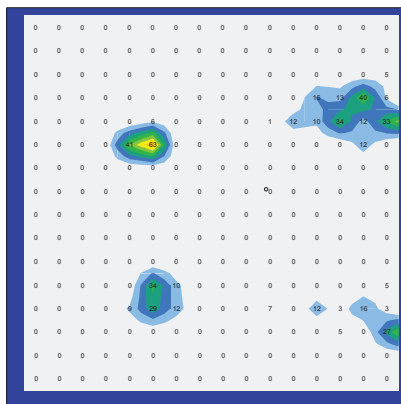

Minimum (mmHg)	0.00
Maximum (mmHg)	63.30
Average (mmHg)	1.86
Variance (mmHg ²)	51.88
Standard deviation (mmHg)	7.20
Coefficient of variation (%)	386.45
Horizontal center (in)	10.99

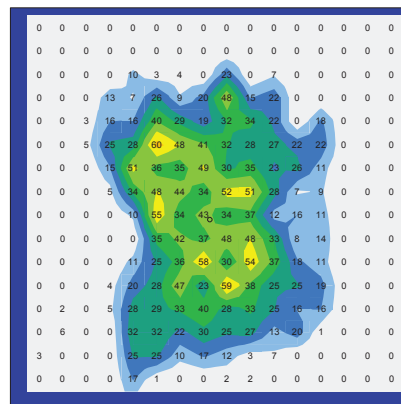

Minimum (mmHg)	0.00
Maximum (mmHg)	59.59
Average (mmHg)	11.98
Variance (mmHg ²)	261.73
Standard deviation (mmHg)	16.18
Coefficient of variation (%)	135.07
Horizontal center (in)	8.26

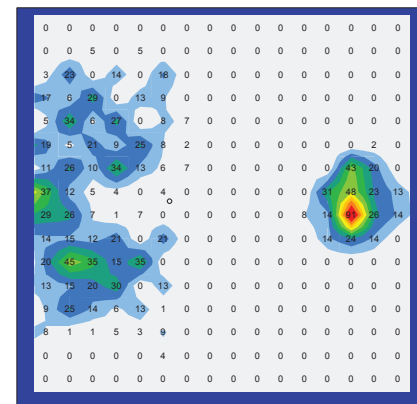

Minimum (mmHg)	0.00
Maximum (mmHg)	91.11
Average (mmHg)	5.33
Variance (mmHg ²)	120.21
Standard deviation (mmHg)	10.96
Coefficient of variation (%)	205.61
Horizontal center (in)	6.13

