



BARIATRIC SCALES

featuring
DIGITAL
Height Rods

PROMED

DETECTO
Model 6857DHR

**DIGITAL SPEED
AND ACCURACY**

Heavy-duty, yet lightweight, the 6857DHR bariatric scale's balanced design and integral wheels make it easy to move. The 6857DHR is engineered specifically for weighing obese and unsteady patients and features a roomy, spacious 2-ft-square platform.



DETECTO DIGITAL BARIATRIC SCALES

-- WITH DIGITAL HEIGHT RODS --



ProMed®

DIGITAL WEIGHT • HEIGHT • BMI

The 6857DHR and 6854DHR ProMed® bariatric scales are designed for optimum performance and display weight, height, and BMI simultaneously (shown below) by

VERSATILE CONNECTIVITY

Detecto's state-of-the-art 758C digital weight indicator includes DB9 serial, USB, or optional Ethernet output for data transfer to a printer or a PC-based EMR (electronic medical record) software program. With the time/date feature, the scale can send the weight, time, date, BMI (Body Mass Index), and height directly to a printer or a PC for record keeping.

simply raising the digital height rod while a patient is being weighed on the scale. There are no keys to push, no operational sequences to memorize, it's fully automatic. These bariatric scales are EMR-ready with their USB and serial ports. The ability to print or save weight data with these scales is a real time-saving device for clinics, doctor's offices, long-term care facilities, and other medical institutions.



- ▶ Digital measurement speed and accuracy
- ▶ Integral wheels for quick and easy transport
- ▶ Stainless steel handrails for assurance
- ▶ Height rod measures up to 6'7" / 200 cm

DETECTO DIGITAL BARIATRIC SCALE SPECIFICATIONS

	ProMed® 6857DHR	ProMed® 6854DHR
CAPACITIES:	800 lb x 0.2 lb / 360 kg x 0.1 kg	600 lb x 0.2 lb / 270 kg x 0.1 kg
PLATFORM SIZE:	24" x 24" / 61 cm x 61 cm	18" x 14" / 46 cm x 36 cm
SHIPPING WEIGHT:	85 lb / 39 kg	74 lb / 34 kg
DIGITAL HEIGHT ROD MEASUREMENTS:	43" / 110 cm - 79" / 200 cm	same
POWER:	6 "C" size alkaline batteries (not included) or optional AC adapter	same
DISPLAY:	Liquid Crystal, 0.5"/13 mm high weight readouts	same
CONNECTIVITY PORTS:	Serial, USB, and optional Ethernet	same




DETECTO
Model 6857DHR



DETECTO
Model 6854DHR

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

Sold By:

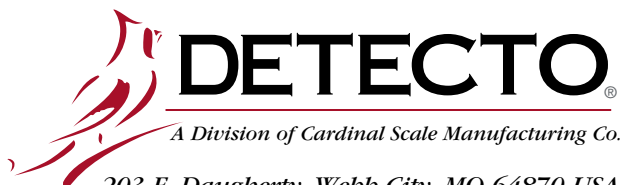

DETECTO®

A Division of Cardinal Scale Manufacturing Co.

203 E. Daugherty, Webb City, MO 64870 USA

Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001

www.Detecto.com