



# The LEGACY

Dolomite developed the Legacy to accommodate a wide variety of users with the most standard features for the best value. The design incorporates these features: flip-up seat with security latch, folding frame, adjustable handles, and wide wheels. With choice of color and an extensive accessory line, the Legacy looks smart strolling down the sidewalk.

## *Frame*

The strongest frame on the market, Dolomite welds the steel frame for greater durability, no screws to loosen at joints. Legacy is available in two colors: Metallic Blue or Candy Apple Red. There are three frame sizes to accommodate individuals of all sizes. Dolomite frames are easily folded for storage or transportation. The handles are height-adjustable in 4 increments. Frame is warranted for 2 years against manufacturer's defects.

## *Braking*

Legacy comes with the patented Dolomite Automatic Braking System with adjustable tension for varied grip strength. Brake handles have wide flange and large loops for comfortable hand positioning. When parked and using the seat, push down the handle to set the parking brakes. The brake mechanism is located within fork for protection from damage. Accessory slow down brake provides adjustable tension for hilly terrain or to even out direction of travel for stroke patients. Warranty is 5 years on brake handle assembly.

## *Standard Features*

Legacy supports 350 pounds, which is 75 pounds more than the majority of the other rolling walkers. Other standard features are:


- Removable basket with carry handle
- Flip-up seat for greater reach or gait training
- Seat handle releases to fold walker
- Magnetic seat latch secures against accidental closure
- Foam grips for comfortable and secure grip
- Flexible and adjustable back strap
- 8" indoor/outdoor tires
- Curb climber on rear wheel assists in lifting the front of walker



## **CLARKE HEALTH CARE PRODUCTS, INC.**

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

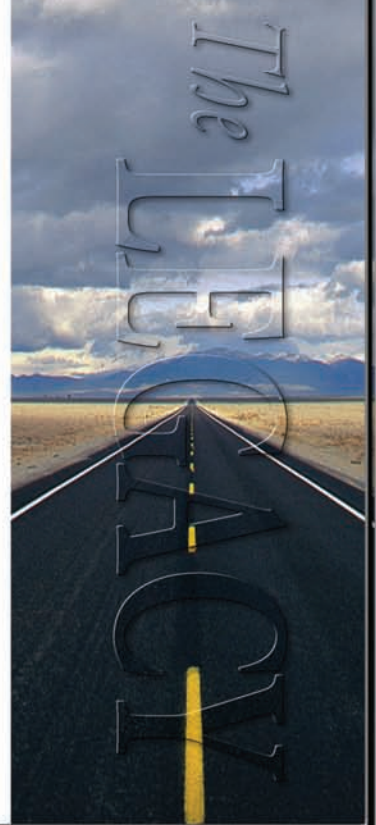



# The LEGACY

## Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



FOOD & BEVERAGE TRAY



CANE HOLDER



TALL HANDLES



UNDER SEAT BAG

Accessories can enhance the use of your walker – the Tray gives you a platform to carry food, the Oxygen Tank Holder will safely support a B, C or D tank. A Cane Holder can secure a cane or crutch to the side of the walker, Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes can control your braking from one handle. Under Seat Bag gives you a place to stash items out of sight. Tall handles add 3" to the handle height.

SPECIFICATIONS	SUPER LOW	LOW	STANDARD
Dist. Between Rear Wheels	19"	19"	19"
Minimum Handle Height	28"	31"	33"
Maximum Handle Height	31"	34"	37"
Minimum & Maximum Tall Handle Heights	31 - 35"	33 - 37"	35 - 39"
Seat Height	19"	22"	24"
Overall Width	23"	23"	23"
Overall Length	23"	26"	26"
Width Between Handles	17"	17"	17"
Folded Dimensions	36" Length x 23" Width x 12" Height		
Walker Weight (lbs.)	15	15	15
Weight Capacity (lbs.)	350	350	350
Product Number: Blue	D12054-37	D12052-37	D12050-37
Product Number: Red	D12054-38	D12052-38	D12050-38

### CLARKE HEALTH CARE PRODUCTS, INC.

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

