


Sleep freely on pillows of air









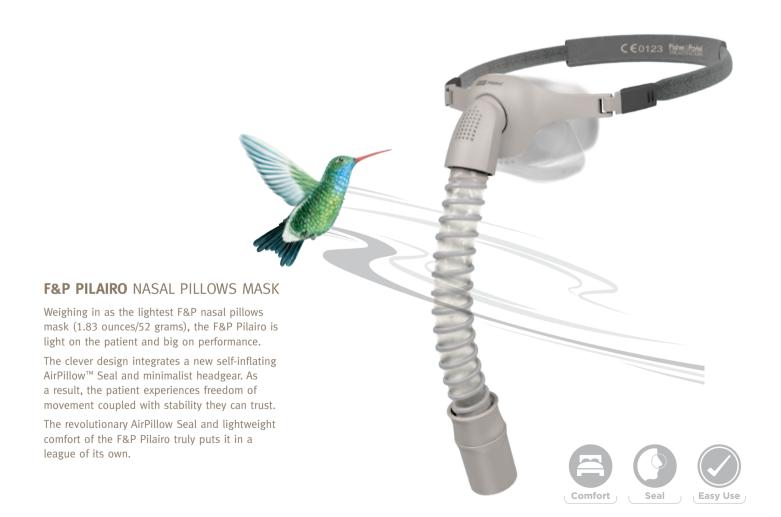



Every detail of the F&P Pilairo has been carefully researched, developed and engineered to improve patient SleepLife*. The F&P Pilairo represents minimalist, highly effective nasal pillows design at its best and sets new benchmarks for comfort, seal and easy use.





small, fast and clever.

Feel the new AirPillow Seal.

- The micro-fine contoured silicone used in the F&P Pilairo's AirPillow Seal is (0.04 inches/0.12 mm) at its thinnest point and only (0.12 inches/3.0 mm) at its thickest.
- Soft, comfortable and highly effective, this fine seal is ergonomically
 designed to self-inflate with (Continuous Positive Airway Pressure) CPAP flow,
 to gently envelop the nose.

Experience total comfort on pillows of air.

- The F&P Pilairo makes the need for complicated adjustments and complex headgear redundant.
- Its minimalist design is stable, low profile and unobtrusive.
- No more messy straps, adjustments or facial indents.

Experience easier therapy.

- One convenient size.
- Only three simple parts to clean.
- Quick and easy to fit.









A little bird told us to make it light,

The design of the F&P Pilairo was inspired by the aerodynamic flight of the world's lightest little bird – the hummingbird. The hummingbird can rotate its fine wings in a circle; it is the only bird able to fly forwards, backwards, up, down and sideways and can hover in midair. This drove the design of a mask seal which self-inflates around the nares in a 'hover-like way' while providing total freedom of movement.





1. AirPillow Seal

- The AirPillow Seal is 3D contoured with fine Variable Thickness Silicone (VTS™).
- It uses CPAP pressure to self-inflate to create an effective yet gentle seal into and around the nose.
 The AirPillow Seal can move and flex in any direction without compromising the seal.

2. Stabilizers

 The stabilizers provide the support frame for the AirPillow Seal and assist with stability.

3. Stretchwise™ Headgear

- Made from a soft latex-free elasticized thread with intelligent-stretch properties, the headgear automatically customizes to the wearer.
- · No manual adjustments are necessary.
- One convenient size.

4. Flexible Tube

- The delivery tube is lightweight, very flexible and designed to minimize pressure or pull on the face.
- The Ball and Socket Elbow provides greater flexibility for tube positioning.

5. Quiet Air Diffuser

 Efficiently designed to ensure the flow is directed away from the patient and their bed partner for a quiet sleep.

FITTING THE F&P PILAIRO IS EASY. FOLLOW THE SIMPLE SEQUENCE BELOW FOR BEST RESULTS:



1. Spread the Wings

Use your fingers to pull apart the pillows.



2. Position the Mask

Place the pillows on your nose until you feel a comfortable seal.



3. Locate Headgear

Pull on the autoadjusting, one-size stretchwise headgear.



4. Test the Wings

Give the wings (pillows) a little wiggle for comfort.



5. Turn on CPAP device

The F&P Pilairo's AirPillow Seal will selfinflate and seemingly hover in and around your nose.

For more information please contact your local Fisher & Paykel Healthcare representative

Manufacturer

Fisher & Paykel Healthcare Ltd 15 Maurice Paykel Place East Tamaki, Auckland 2013

PO Box 14 348 Panmure Auckland 1741 New Zealand

+64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.co.nz Web: www.fphcare.com

Australia

Fisher & Paykel Healthcare Pty Limited 36-40 New Street, PO Box 167 Ringwood, Melbourne Victoria 3134. Australia

Tel: +61 3 9879 5022 Fax: +61 3 9879 5232

Austria

Tel: 0800 29 31 23 Fax: 0800 29 31 22

Benelux

Tel: +31 40 216 3555 Fax: +31 40 216 3554

Tel: +86 20 3205 3486 Fax: +86 20 3205 2132

France

Tel: +33 1 6446 5201 Fax: +33 1 6446 5221

Germany

Tel: +49 7181 98599 0 Fax: +49 7181 98599 66

India

Tel: +91 80 4284 4000 Fax: +91 80 4123 6044

Irish Republic

Tel: 1800 409 011

Italy

Tel: +39 06 7839 2939 Fax: +39 06 7814 7709

Japan

Tel: +81 3 3661 7205 Fax: +81 3 3661 7206

Northern Ireland

Tel: 0800 132 189

Spain

. Tel: +34 902 013 346 Fax: +34 902 013 379

Tel: +46 8 564 76 680 Fax: +46 8 36 63 10

Switzerland

Tel: 0800 83 47 63 Fax: 0800 83 47 54

Taiwan

Tel: +886 2 8751 1739 Fax: +886 2 8751 5625

Turkey

Tel: +90 312 354 34 12 Fax: +90 312 354 31 01

UK EC REP

Fisher & Paykel Healthcare Ltd Unit 16, Cordwallis Park Clivemont Road, Maidenhead Berkshire SL6 7BU, UK

Tel: +44 1628 626 136 Fax: +44 1628 626 146

USA/CANADA

Tel: +1 800 446 3908 or +1 949 453 4000 Fax: +1 949 453 4001

フィッシャー&パイケルヘルスケア株式会社

7103-0025 東京都中央区日本橋茅場町3-8-10 リベラ茅場町5F FAX. 03-3661-7206 ご注文書専用 FAX. 0800-100-7207 TEL. 03-3661-7205 (代表)

- · Complete mask with headgear assembled and ready for use
- · One convenient size
- User-friendly tear off Fitting Guide

Reorder Code

400420 F&P Pilairo

